

NITgold Cit 30nm

#510010 & #510011

STORE AT 4°C away from light. DO NOT FREEZE

Description

The 30nm gold nanoparticles are uniform quasi-spherical gold nanoparticles with a mean diameter of thirty nanometers capped with citrate. These nanoparticles present unique physic-chemical properties (i.e. optical, catalytic, electrochemical activity, etc.). Citrate gold nanoparticles are excellent candidates to be functionalized since a wide range of molecules can be tethered onto the AuNP surface by means of a thiol (SH) group displacing adsorbed citrate anions. Once functionalized, they can be employed as platforms for many applications such as target-specific drug delivery, sensors, lateral flow tests, imaging probes for dark-field microscopy, flow cytometry, cancer photothermal therapy, catalysis, and optoelectronic. They can be also conjugated with distinct molecules and/or other micro/nanostructures to generate new biosensors, therapeutics, etc. The optical and electronic properties of gold nanoparticles are tuneable by changing their size, shape, surface chemistry or even their aggregation state.

Technical Specifications

Particle Surface: Citrate capping

Peak SPR wavelength: 526 nm

Average Diameter:¹ 29.3 ± 3.5 nm

Hydrodynamic Diameter (DLS): 30.1 ± 3.2 nm

Molar Concentration:² 0.28 nM

O.D.: 1

Particle Concentration: 1.7×10¹¹ particles/mL


Z-Potential: -36.2 mV

Solvent: Milli-Q Water


pH of Solution: 5.0

¹ Particle diameter determined by Surface Plasmon Resonance from the spectrum. Haiss et al. Anal. Chem. 2007, 79, 4215-4221.


² Determined according to a $\epsilon = 1.74 \times 10^9 \text{ M}^{-1} \cdot \text{cm}^{-1}$. Werts et al. Analyst 2013,138, 583-592.

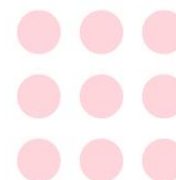

UV/visible absorbance spectrum

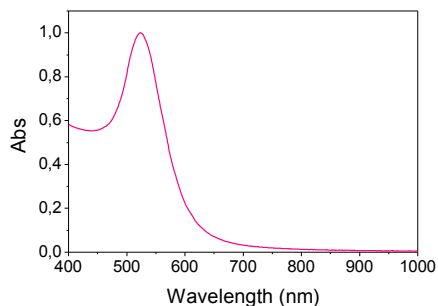

Size Distribution (DLS)

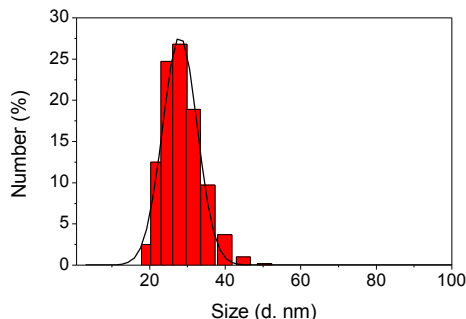

Suggested Application(s)

- Biosensing
- Catalysis
- Drug delivery
- Colorimetric probes
- Lateral flow immunoassays

Ordering Information

Product Name	Nº Nanoparticles/mL	O.D.	ϵ (M ⁻¹ cm ⁻¹)	Quantity	Catalogue No.
NITgold Cit 30nm	1,71E+11	10D	3,52E+09	25mL	51001017V
NITgold Cit 30nm	1,71E+11	10D	3,52E+09	100mL	510010160
NITgold Cit 30nm	1,71E+11	10D	3,52E+09	250mL	51001016C
NITgold Cit 30nm	1,71E+11	10D	3,52E+09	500mL	51001016F
NITgold Cit 30nm	5,13E+11	30D	3,52E+09	25mL	51001117V
NITgold Cit 30nm	5,13E+11	30D	3,52E+09	100mL	510011160
NITgold Cit 30nm	5,13E+11	30D	3,52E+09	250mL	51001116C

Product disclaimer

This nitparticles® product is to be used for research purposes only. Unless stated in the documentation of on an individual product label, catalog or other information provided to the buyer, IT IS FORBIDDEN TO USE IT for different purposes, including but not limited to them: in vitro diagnostic, use in food, pharmaceutical purposes, medical purposes, or use in cosmetic products, neither for use in humans nor animals, nor for any commercial purposes. Please refer to www.nitparticles.com for the Material Safety Data Sheet of the product.

